

REINFORCED FLEX SLEEVE RFX

Gel Suction Suspension Sleeve

ALPS® Reinforced Flex prosthetic suspension sleeve delivers the same qualities as our popular SFX but with an extra layer of durability. While the OptiGel and grey knitted fabric combine to supply maximum comfort and high performance, the RFX also features an added layer of internal fabric to offer increased durability around the trim line of the user's socket. At the same time, the sleeve provides a 30° pre-flexion to ease bending and reduce bunching behind the knee, making the RFX an ideal choice for active users.

 Available in 3mm uniform thickness.



TOP BENEFITS

- ✓ 30° pre-flexion to ease bending and reduce bunching behind the knee
- ✓ Provides a secure suspension
- ✓ OptiGel offers a great degree of control
- ✓ Available in one length that accommodates most users


Recommended for K3-K4 activity levels.



TOP FEATURES


- ✓ Extra layer of internal fabric provides a more durable option for high-activity users
- ✓ Grey knitted fabric combined with OptiGel delivers maximum comfort and high performance

CROSS SECTION DIAGRAM

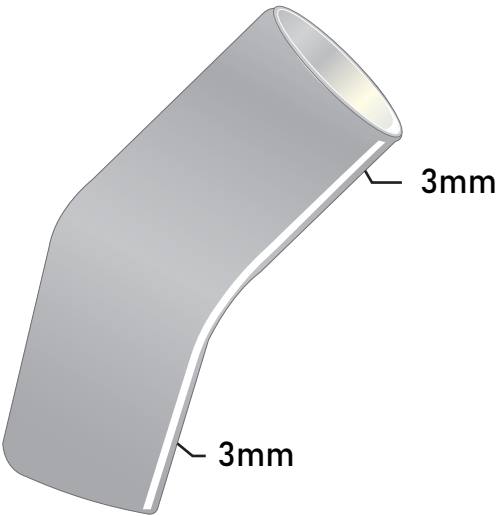

HOW TO SIZE

Measurements taken 10cm above the patella

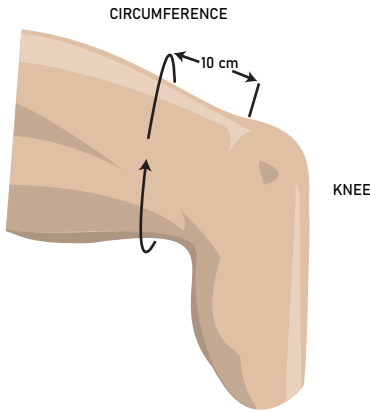

SIZING CHART

MEASURED SIZE	3.0 MM UNIFORM
34-37 cm	RFX20-3G
36-39 cm	RFX23-3G
38-42 cm	RFX26-3G
41-46 cm	RFX29-3G
45-50 cm	RFX32-3G

Blatchford distribute these products on behalf of ALPS® who have supplied all technical content for this publication.
Please refer to the product IFU for additional information.